

Office européen des brevets

(11) Publication number:

0 184 056  
A3

(12)

EUROPEAN PATENT APPLICATION

(21) Application number: 85114561.5

(51) Int. Cl.: C 07 H 21/00

(22) Date of filing: 16.11.85

INF 0201 0511

(30) Priority: 27.11.84 US 675386

(71) Applicant: Molecular Diagnostics, Inc., 400 Morgan Lane, West Haven, CT.06516 (US)

(43) Date of publication of application: 11.06.86  
Bulletin 86/24

(72) Inventor: Dattagupta, Nanibhushan, Dr., 470 Prospect Street, New Haven, CT 06511 (US)  
Inventor: Rae, Peter, Dr., 71 Ingram Street, Hamden, CT 06517 (US)  
Inventor: Crothers, Donald, Dr., 8 Surrey Lane, Northford, CT 06472 (US)  
Inventor: Barnett, Thomas, Dr., 27 Jeffrey Rd, East Haven, CT 06512 (US)

(84) Designated Contracting States: AT BE CH DE FR GB IT LI NL SE

(83) Date of deferred publication of search report: 15.04.87 Bulletin 87/16

(74) Representative: Jesse, Ralf-Rüdiger, Dr. et al, Bayer AG Konzernverwaltung RP Patentabteilung, D-5090 Leverkusen 1 Bayerwerk (DE)

(54) Large scale production of DNA probes.

≠ 11

(57) A process for production of a single strand of a nucleic acid comprising covalently linking to a solid substrate a polynucleotide complementary to the desired strand, hybridizing said polynucleotide with an oligonucleotide, extending the oligonucleotide in direction away from said substrate, denaturing the hybridized polynucleotide and extended oligonucleotide, thereby to free the extended oligonucleotide from the solid substrate, and separating the extended oligonucleotide. The product can be used for making analytical probes.

*Primer elongation*

FOI290 906E960

EP 0 184 056 A3

19



Europäisches Patentamt  
European Patent Office  
Office européen des brevets

11

Publication number:

0 184 056  
A2

12

## EUROPEAN PATENT APPLICATION

21

Application number: 85114561.5

51

Int. Cl. 4: C 07 H 21/00


22

Date of filing: 16.11.85

No	références, formules, pages à photocopier, etc	No	classement
	DOC	1	C07H21/00 C4
		2	C12Q1/68



11.01.10.93


 Europäisches Patentamt  
 European Patent Office  
 Office européen des brevets

(11) Publication number:

0 184 056  
A2

(12)

# EUROPEAN PATENT APPLICATION

(21) Application number: 85114561.5

(51) Int. Cl.<sup>4</sup>: C 07 H 21/00

(22) Date of filing: 16.11.85

(30) Priority: 27.11.84 US 675386

(43) Date of publication of application:  
11.06.86 Bulletin 86/24

(84) Designated Contracting States:  
AT BE CH DE FR GB IT LI NL SE

(71) Applicant: Molecular Diagnostics, Inc.  
400 Morgan Lane  
West Haven, CT.06516(US)

(72) Inventor: Dattagupta, Nanibhushan, Dr.  
470 Prospect Street  
New Haven, CT 06511(US)

(72) Inventor: Rae, Peter, Dr.  
71 Ingram Street  
Hamden, CT 06517(US)

(72) Inventor: Crothers, Donald, Dr.  
8 Surrey Lane  
Northford, CT 06472(US)

(72) Inventor: Barnett, Thomas, Dr.  
27 Jeffrey Rd  
East Haven, CT 06512(US)

(74) Representative: Jesse, Ralf-Rüdiger, Dr. et al,  
Bayer AG Konzernverwaltung RP Patentabteilung  
D-5090 Leverkusen 1 Bayerwerk(DE)

(54) Large scale production of DNA probes.

(57) A process for production of a single strand of a nucleic acid comprising a covalently linking to a solid substrate a polynucleotide complementary to the desired strand, hybridizing said polynucleotide with an oligonucleotide, extending the oligonucleotide in direction away from said substrate, denaturing the hybridized polynucleotide and extended oligonucleotide, thereby to free the extended oligonucleotide from the solid substrate, and separating the extended oligonucleotide. The product can be used for making analytical probes.

11.01.10.93

EP 0 184 056 A2

This invention relates to a method for the large scale production of a specific nucleic acid sequence:

Large scale production of a nucleic acid sequence is usually done by cloning the particular sequence in a specific vector. The sequence to be cloned is isolated, identified and then coupled covalently to a single or double-stranded vector. The vectors with the extra DNA are separated from the host cell and, depending on the requirements, the cloned piece of DNA has to be restricted and separated from the rest of the DNA. If one requires single-stranded DNA, either it is cloned in a single-stranded vector or strand separator is necessary. All these techniques involve skilled manipulation of biochemical and biological systems.

Analytical Biochemistry, 140, 95-103(1984) describes a method of producing DNA hybridization probes by non-specifically immobilizing single strand DNA template to cellulose. Although the method is useful, the length distribution of the newly synthesized product DNA is not as uniform as might be desired. It now appears this may be due to multiple attachments of the template DNA to the cellulose.

It is accordingly an object of this invention to synthesize specific nucleic acid sequences on a relatively large scale without the continual need for plasmids, cloning and restriction.

This and other objects and advantages are realized in accordance with the present invention pursuant to which there is provided a process for production of a single strand of a nucleic acid comprising covalently linking to a solid substrate a polynucleotide complementary to the desired strand, hybridizing said polynucleotide with an oligonucleotide, extending the oligonucleotide in direction away from said substrate, denaturing the hybridized polynucleotide and extended oligonucleotide, thereby to free the extended oligonucleotide from the solid substrate, and separating

2

the extended oligonucleotide from the solid support carrying the covalently linked polynucleotide which support is then recycled to an earlier step in the process.

5       The intermediate product comprising a solid substrate, a single-stranded polynucleotide covalently linked to the solid substrate, and an oligonucleotide hybridized to the polynucleotide is also new.

10       The invention will be further described with reference to the accompanying drawing which is a flow sheet of a preferred process.

15       A solid support 10 is covalently coupled to a DNA strand 12 complementary to the desired strand to be synthesized, the DNA strand having its 3'-end adjacent to the solid support. This product, 14, is known.

20       An oligonucleotide 16 corresponding to the 5'-end of the desired strand is then hybridized to strand 12 of product 14. This solid new intermediate 18 is contacted with a solution containing nucleoside triphosphates and the Klenow fragment of DNA polymerase, whereupon the oligonucleotide grows at its 3'-end, following the DNA strand 12 serving as a template, producing the desired strand 20 base-paired to 12. The structure constituting the polynucleotide base-paired to the extended  
25       oligonucleotide is denatured so as to release the strand 20 to the solution and then the strand-containing solution is separated from the product 14, as by centrifugation, the product 14 then being recycled to an earlier step in the process.

30       The desired strand 20 can be recovered from its solution in any known manner.

35       Describing the invention in greater detail, the solid support can be cellulose, Sephadex, or Sepharose, a piece of paper, nylon or anything which can be used to react with an amine or aldehyde or similar residue.

      The nucleic acid is linked to the solid support in a known manner, e.g., the latter providing an aldehyde or

amine reactive group which can react with an amine or aldehyde group provided on the nucleic acid to form a Schiff's base, followed by its reduction as with a borohydride to form an aminomethylene joinder between support and nucleic acid. This is described more fully in European Patent Application No. 85 101 407.6. The link can also be a phosphate ester, etc.

After the gene probe or the gene to be copied has been immobilized, it is primed with a suitable piece of oligonucleotide such that the 5' end of the oligonucleotide will hybridize to the 3' end of the immobilized DNA. Since the 3' end of the immobilized DNA is closer to the solid support, the 5' end of the oligonucleotide will also be closer to the solid support and the 3' hydroxyl residue of the oligonucleotide primer will be available for primer extension by several DNA synthesis enzymes, for example, DNA polymerase, Klenow fragment of DNA polymerase or reverse transcriptase. After the primer has been hybridized to the immobilized DNA, it can be extended by using those enzymes and suitable nucleoside triphosphate substrates. After the enzyme reaction is over, and the enzyme is washed away, one can separate the newly synthesized DNA strand by heating the solution or by adding some denaturing agent, for example, sodium hydroxide solution, formamide or dimethyl sulfoxide.

The solid support carrying immobilized nucleic acid is separated from the solution containing the desired strand as by filtration and/or centrifugation and the solid support can be recycled, if desired, after washing.

The supernatant containing the desired strand is subjected to dialysis, alcohol precipitation or the like, to obtain pure strand which can then be used in known manner. One such use is for making probes for diagnostic

tests as described more fully in European Patent Application No. 84 107 248.1.

The invention will be further described in the following illustrative examples wherein all parts are by weight unless otherwise expressed or apparent from the context.

#### Example 1

Site Specific Immobilization of the DNA for the purpose of large scale production of the DNA probes for sickle cell anemia diagnosis.

A segment of  $\beta$ -globin gene containing 737 nucleotide residues has been cloned in a vector M13mp7. Since the cloning sites on M13mp7 are in an inverted repeat, they exist in a double-stranded form even in single stranded M13 DNA. This double-stranded site can be digested with proper restriction enzyme. M13mp7 vector with  $\sim$ 740 residues of  $\beta$ -globin inserted is first digested with the restriction enzyme which is either BamH1 or EcoR1. The strands are then separated by gel electrophoresis. The fragment containing the globin insert is immobilized onto Sephadex G25 by the following procedure:

First the Sephadex particle is activated with CNBr as has been described in Example 2. The 737 globin insert-containing fragment is tailed with terminal deoxynucleotidyl transferase using ATP as a substrate (Examples 3 and 4). ATP, being a ribonucleoside triphosphate, contains a sugar residue which can be oxidized to produce a dialdehyde. These dialdehyde residues can then be reacted with the amide on the solid support to form a Schiff base which can be reduced to produce secondary amine linkage between the solid support and the DNA.

After the probe is immobilized, a primer containing 19 nucleotides is hybridized with the immobilized DNA and then the primer extension is performed by using either

Klenow fragment of DNA polymerase I or reverse transcriptase as will now be described.

Example 2

5        Preparation of the solid support containing Amine  
       (NH<sub>2</sub>) groups.

      Sephadex G25 powder (0.5 g) is mixed with 5 ml  
       distilled water at room temperature for approximately 30  
       minutes. The product is washed with distilled water and  
 10        resuspended in 20 ml ice cold distilled water. The  
       suspension is adjusted to pH 10.5 - 11 with a 5 M  
       solution of sodium hydroxide. The suspension is cooled  
       in an ice-bath (4°C) and 1.0 g solid cyanogen bromide is  
       added to the suspension. The suspension is retained in  
 15        the ice-bath and continuously stirred for 30 minutes.  
       During this period, the pH is continuously monitored and  
       maintained at 10.5 - 11 with sodium hydroxide (5 M). The  
       suspension is washed with ice cold distilled water and  
       treated with approximately 5 ml of hexamethylene diamine  
 20        (1 M solution in H<sub>2</sub>O, adjusted to pH 9 - 10 with HCl).  
       The Sephadex amine product is washed with water or other  
       suitable buffer, depending on the future use of the  
       product.

25        Example 3

Terminal transferase addition of ribonucleotides.

      A ribonucleotide or deoxyribonucleotide can be  
       coupled to the 3'-hydroxyl terminus of a DNA molecule or  
       oligonucleotide by the enzymatic action of terminal  
 30        deoxynucleotidyl transferase. Twenty (20) µl DNA 737 bp  
       approximately  $5 \times 10^{-13}$  moles are dissolved in 20 µl  
       potassium cacodylate (1 M, pH 7.2). This solution is  
       incubated with 40 µl distilled water, 5-20 µl dithio  
       threitol (DTT) (2 mM), 1-2 µl ribonucleoside triphosphate  
 35        (10 mM) and 50 µl magnesium chloride (20 mM) at 37°C for  
       5 minutes. The reaction mixture is subsequently cooled  
       in an ice-bath for 10 minutes. Terminal deoxynucleotidyl.

053906 0640



00:01:19 M

0184056

transferase (14 - 20 units)<sup>6</sup> is added and the mixture is incubated at 15° for 24 hours. The modified nucleic acid can be isolated from non-coupled nucleotides and protein by phenol extraction and alcohol precipitation.

5

Example 4Coupling of end labeled nucleic acids to a solid support.a) Coupling via Schiff's base formation.

10 It is known that cis-diols can be oxidized by periodate to form dialdehydes. The dialdehydes can form Schiff's base with a primary amine group via the addition of  $-NH_2$ . The Schiff's base can be reduced with sodium borohydride to form a secondary amine.

15 In a similar manner, (DNA containing a ribonucleotide at the 3'-hydroxyl terminus (Example 3) is dissolved in sodium acetate buffer, 0.1 M, pH 5 at a concentration of 1 mM. Twenty (20)  $\mu$ l of sodium metaperiodate (100 mM) is added to 1 ml of the nucleic acid solution. The reaction is allowed to proceed for 40 minutes at room temperature (25 °C). Following the reaction, the pH is adjusted to 8 with sodium hydroxide solution. The solution is added to a Sephadex-amine (Example 2), prepared as described above, suspended in a  
25 suitable buffer such as 0.1 M sodium acetate pH 8. The reaction is allowed to proceed for 30 minutes at room temperature (25 °C) resulting in the formation of a suspension of Schiff's base. The Schiff's base is reduced by the addition of sodium borohydride. The  
30 reduction is carried out in four steps: Approximately 0.15 ml of freshly prepared sodium borohydride solution (200 mM) is added and the reaction is allowed to proceed for 30 minutes. Approximately 0.15 ml of sodium borohydride solution is again added to the reaction  
35 mixture and the reaction is continued for 60 minutes. Another 0.15 ml of the sodium borohydride is subsequently added to the reaction mixture. After 90 minutes, another

MD 2110.93 0184055  
aliquot of 0.15 ml of the sodium borohydride solution is added to complete the reaction. The suspension is centrifuged for 30 minutes at 2000 g. The supernatant is decanted and the wash is repeated 3 times. The pellet is resuspended in 5 ml of the hybridization buffer as in Example 5.

#### Example 5

Hybridization of the oligonucleotide 19A' with the immobilized 737 insert as in Example 4:

The 737 insert immobilized support is suspended in a buffer containing 50 mM tris and 1 mM dithiothreitol. Oligonucleotide 19A', Proc. Natl. Acad. Sci. USA 80, 278 (1983) is dissolved in the same buffer and is added to the suspension. The whole mixture is heated to 65°C and then brought down slowly to 50°C. The mixture is left for 15 minutes after which time it is cooled to room temperature. Unreacted, unhybridized oligonucleotide is removed by washing at 0°C with the same buffer. The hybrid containing the immobilized support DNA is washed with 15 mM tris buffer, pH 7.2.

#### Example 6

The preparation of the hybridization probe using reverse transcriptase of the Klenow fragment of DNA polymerase 1.

The hybrid immobilized onto the solid support produced in Example 2 is taken in a buffer containing 50 mM tris, pH 7.2, 10 mM magnesium sulfate, 1 mM dithiothreitol, 50 micrograms per ml BSA and 200 nM each of 1-deoxy ATP, deoxy CTP, deoxy GTP, deoxy TTP and, e.g., 5 units reverse transcriptase or Klenow fragment of DNA polymerase. The whole mixture is incubated at room temperature for 30 minutes and the reaction is stopped by centrifuging the supernatant and removing the solid particles. The solid particles are washed with the same buffer and heated to 40°C after adding 80% formamide to

09.01.1974

0184056

remove the newly synthesized probe from the solid support. The support is washed thoroughly with water, then with tris buffer and it is recycled as in Example 2 and 3. The synthesized probe is then analyzed by gel electrophoresis to determine the length.

Example 7

Cloning of a beta-globin gene segment in M13

To prepare a single-stranded segment of the beta-hemoglobin gene that can be separated from the phage DNA in which it is propagated, an approximately 740 base pair Alu I segment of the globin gene was inserted into the linker region of the double stranded replicative form of M13mp7.

The linker region that has been incorporated into the phage DNA is not only rich in restriction enzyme cleavage sites, but is a reverse repeat, and has the priority of forming a duplex in otherwise single-stranded phage DNA. Restriction enzymes require double-stranded DNA as a substrate; thus the mp7 phage DNA can be cleaved by the enzyme. Cloning of a foreign DNA segment into the Hinc II site, for example, allows digestion of the single-stranded phage DNA by such enzymes as EcoRI or Bam HI, separating the recombinant phage DNA into vector and insert strands that are readily separable from one another.

The 740 bp Alu I segment of the beta-globin gene comprises the beginning of the gene and a portion of its upstream flank. The segment was initially cloned by converting the Alu I ends to EcoRI ends, and inserting the segment into pBR322 at its EcoRI site. Recloning the EcoRI segment into M13mp7 in a useful location involved the filling-in of EcoRI ends of the globin segment so that it could be blunt-end ligated into mp7 DNA digested with Hinc II. Insertion at the Hinc II site assured that the

insert could be effectively removed from single-stranded recombinant phage DNA by EcoRI or Bam HI digestion.

Establishment of the orientation of the insert in clones, and determination of whether the insert strand in phage was coding or noncoding, was accomplished by restriction enzyme digestion, and hybridization of strand-specific synthetic probes to phage DNA.

#### Example 8

The newly synthesized probe of Example 6 hybridizes with the 201 bp fragment of the Dde I digested nucleic acids containing  $\beta$ -globin gene. This probe is used for the detection of the digestion pattern in a sickle and normal  $\beta$ -globin genomic DNA following the procedure described in U.S. Patent No. 4,395,486.

It will be understood that the specification and examples are illustrative but not limitative of the present invention and that other embodiments within the spirit and scope of the invention will suggest themselves to those skilled in the art.

1. A structure for producing a specific nucleic acid strand, comprising a solid substrate, a single-stranded polynucleotide covalently linked at one end to the solid substrate, and an oligonucleotide hybridized to the polynucleotide.

2. A structure according to claim 1, wherein the oligonucleotide comprises at least five nucleotide residues.

3. A structure according to any of claims 1 and 2, wherein the 3'-end of the polynucleotide is linked to the solid substrate.

4. A structure according to any of claims 1 to 3, wherein the link contains a  $\text{-}\overset{|}{\underset{|}{\text{N}}}\text{-}\overset{|}{\underset{|}{\text{C}}}\text{-}$  bond.

5. A structure for producing a specific nucleic acid strand, comprising a solid substrate, a single-stranded first polynucleotide covalently linked at one end to the solid substrate, and a second polynucleotide hybridized to the first polynucleotide, one end of one polynucleotide being complementary to the other end of the second polynucleotide.

6. A structure according to claim 5, wherein the 3'-end of the first polynucleotide is linked to the solid substrate, the link containing a  $\text{-}\overset{|}{\underset{|}{\text{N}}}\text{-}\overset{|}{\underset{|}{\text{C}}}\text{-}$  bond.

7. A process for production of a single strand of a nucleic acid comprising covalently linking to a solid substrate a polynucleotide complementary to the desired strand, hybridizing said polynucleotide with an oligonucleotide, extending the oligonucleotide in

direction away from said substrate, denaturing the hybridized polynucleotide and extended oligonucleotide, thereby to free the extended oligonucleotide from the solid substrate, and separating the extended oligonucleotide.

5

8. A process according to claim 7, including the further step of recycling the solid support carrying the covalently linked polynucleotide.

10

15

20

25

30


35

0953906 062401

M 21.10.93

0184056

1/1

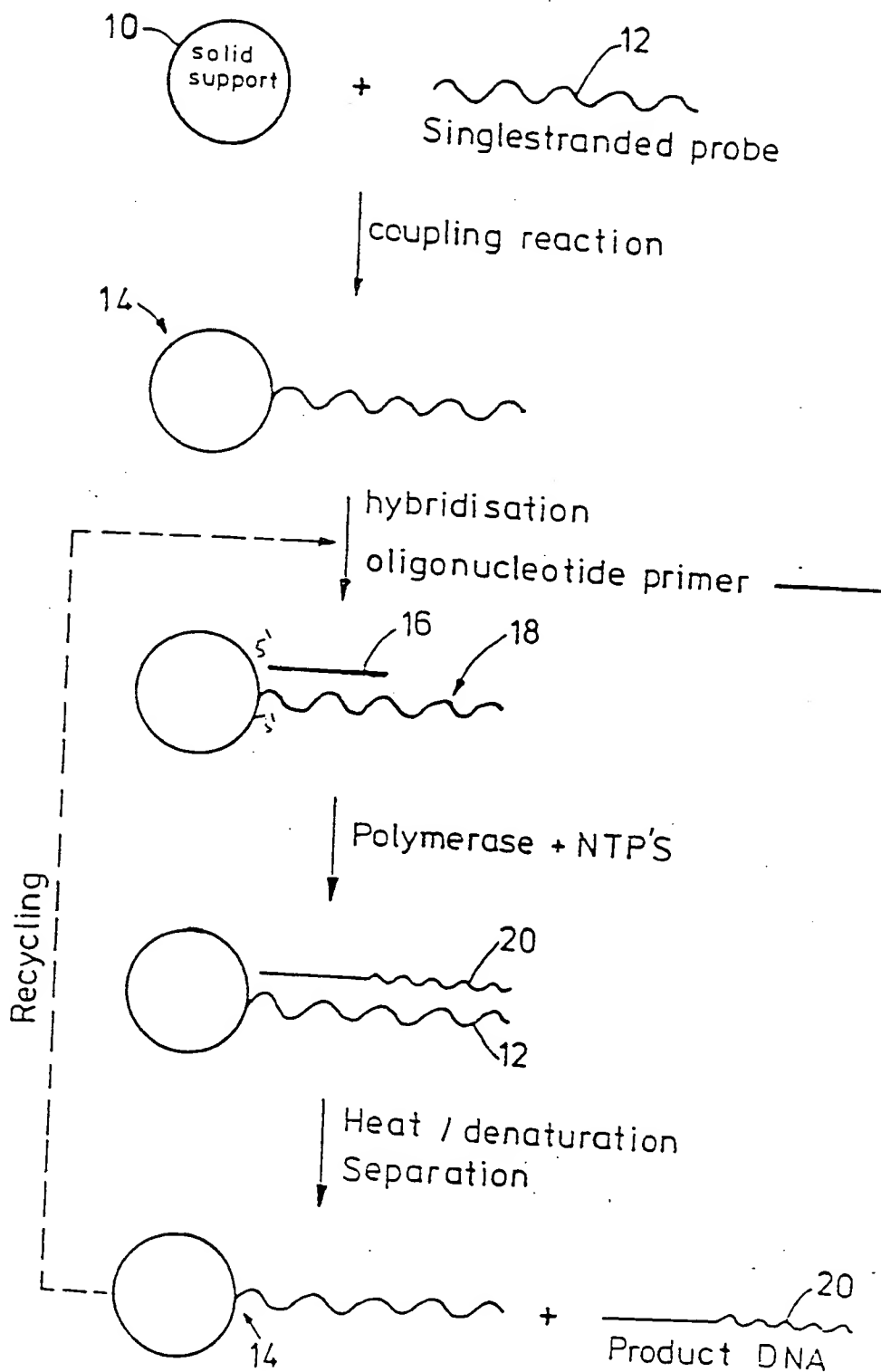

09.01.1986



European Patent  
Office

EUROPEAN SEARCH REPORT

0184056

Application number

EP 85 11 4561

DOCUMENTS CONSIDERED TO BE RELEVANT

Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int. Cl.+)
	No relevant documents have been disclosed  -----		C 07 H 21/00
			TECHNICAL FIELDS SEARCHED (Int. Cl.+)
			C 07 H 21/00 C 12 N 15/00
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 16-12-1986	Examiner VAN BIJLEN H.

CATEGORY OF CITED DOCUMENTS

- X : particularly relevant if taken alone  
Y : particularly relevant if combined with another document of the same category  
A : technological background  
D : non-written disclosure

- T : theory or principle underlying the invention  
E : earlier patent document, but published on, or after the filing date  
O : document cited in the application  
L : document cited for other reasons  
& : member of the same patent family, corresponding

09.01.1986



Europäisches Patentamt,  
Erhardtstraße 27,  
W-8000 München 2

European Patent Office,  
Erhardtstrasse 27,  
W-8000 Munich 2

Office européen des brevets,  
Erhardtstr. 27,  
W-8000 Munich 2

An alle Patentanmelder, Patent-  
inhaber und ihre Vertreter

To all applicants, patent  
proprietors and representatives

A tous les demandeurs, les titu-  
laires de brevet et leurs manda-  
taires

## W A R N U N G

## W A R N I N G

## AVERTISSEMENT

vor irreführenden Anpreisungen  
und Zahlungsaufforderungen  
von Firmen, die angeblich euro-  
päische Patente registrieren

Beware of approaches and  
requests for payment from firms  
purporting to register European  
patents

concernant les offres fraudu-  
leuses et les invitations à payer  
émanant de sociétés qui pré-  
tendent procéder à l'enregistre-  
ment de brevets européens

Zahlreiche Patentanmelder und Pa-  
tentinhaber haben das Europäi-  
sche Patentamt darauf aufmerk-  
sam gemacht, daß Unternehmen  
und Einzelpersonen - unter Ver-  
wendung amtlicher oder amtlich  
klingender Bezeichnungen - zur  
gesetzlich nicht wirksamen Eintra-  
gung von Patenten in nichtamtli-  
chen Registern oder Publikationen  
einladen. Hierfür werden in der Re-  
gel hohe Geldzahlungen verlangt,  
ohne daß die Patentanmelder bzw.  
-inhaber dafür eine adäquate Ge-  
genleistung erhalten.

Numerous applicants and patent  
proprietors have informed the  
European Patent Office that firms  
and individuals using  
official-sounding names have been  
inviting them to register patents in  
unofficial registers or publications.  
High fees are generally charged,  
for which the patent applicant or  
proprietor receives no real service.

De nombreux demandeurs et titu-  
laires de brevet ont attiré l'atten-  
tion de l'Office européen des bre-  
vets sur le fait que des sociétés et  
des particuliers, utilisant des dé-  
nominations officielles ou à con-  
sonance officielle, offrent d'effec-  
tuer l'inscription - qui est dénuée  
d'effet légal -, de brevets dans des  
registres ou des publications non  
officielles et réclament générale-  
ment des sommes importantes,  
sans contrepartie adéquate pour le  
demandeur ou le titulaire du bre-  
vet.

Das Europäische Patentamt weist  
darauf hin, daß eine gesetzlich  
wirksame Veröffentlichung bzw.  
Eintragung von europäischen Pa-  
tentanmeldungen und Patenten  
nur über das Europäische Patent-  
amt erreicht werden kann und daß  
die Eintragung in das Europäische  
Patentregister von Amts wegen er-  
folgt.

The European Patent Office would  
point out that it alone performs  
legally effective publication and  
registration of European patent  
applications or patents, and that it  
does so automatically.

L'Office européen des brevets sig-  
nale que la publication et l'inscrip-  
tion de demandes de brevet europ-  
éen et de brevets européens ne  
peuvent produire d'effet légal que  
si elles sont effectuées par l'inter-  
médiaire de l'Office européen des  
brevets, et que l'inscription au Re-  
gistre européen des brevets a lieu  
d'office.

Das Europäische Patentamt  
warnt eindringlich vor den irre-  
führenden Machenschaften sol-  
cher Unternehmen und Einzel-  
personen. Insbesondere wird  
darauf aufmerksam gemacht,  
daß die angebotenen Leistun-  
gen keinerlei Rechtswirkung  
entfalten: Eine Zahlungsver-  
pflichtung besteht daher nicht.

The European Patent Office  
warns all concerned to guard  
against the subterfuges of such  
firms and individuals. You  
should also be aware that  
because their services have no  
legal effect whatsoever, there is  
no obligation to pay.

L'Office met instamment en  
garde les demandeurs et leurs  
mandataires contre les agisse-  
ments frauduleux de ces so-  
ciétés et de ces particuliers. Il  
est notamment rappelé que les  
prestations offertes ne produi-  
sent aucune espèce d'effet juri-  
dique, et qu'il n'existe donc au-  
cune obligation de paiement.

Stellen sich für Sie weitere Fragen,  
wenden Sie sich bitte an ihren Ver-  
treter oder an das Europäische Pa-  
tentamt (Tel.: München 089/  
2399 4512).

If in doubt, please contact either  
your representative or the  
European Patent Office  
(Tel. Munich 089/2399 4512).

Pour toute autre question à ce su-  
jet, veuillez vous adresser à votre  
mandataire ou à l'Office européen  
des brevets (tel. Munich 089/2399  
4512).

0953306-062401