

1 COOLEY GODWARD LLP  
STEPHEN P. SWINTON (106398)  
2 J. CHRISTOPHER JACZKO (149317)  
4365 Executive Drive, Suite 1100  
3 San Diego, California 92121-2128  
Telephone: (858) 550-6000  
4 Facsimile: (858) 453-3555

5 DOUGLAS E. OLSON (38649)  
BROBECK, PHLEGER & HARRISON LLP  
6 12390 El Camino Real  
San Diego, California 92130-2081  
7 Telephone: (858) 720-2500  
Facsimile: (858) 720-2555

8 R. WILLIAM BOWEN, JR. (102178)  
9 GEN-PROBE INCORPORATED  
10210 Genetic Center Drive  
10 San Diego, California 92121-4362  
Telephone: (858) 410-8918  
11 Facsimile: (858) 410-8637

12 Attorneys for Plaintiff  
Gen-Probe Incorporated

13  
14 UNITED STATES DISTRICT COURT  
15 SOUTHERN DISTRICT OF CALIFORNIA

16  
17 GEN-PROBE INCORPORATED,

18 Plaintiff,

19 v.

20 VYSIS, INC.,

21 Defendant.

No. 99-CV-2668H AJB  
JUDGE MARILYN L. HUFF

**DECLARATION OF DR. MATTHEW LONGIARU  
IN SUPPORT OF GEN-PROBE'S MOTION FOR  
PARTIAL SUMMARY JUDGMENT**

Date: May 29, 2001  
Time: 10:30 a.m.  
Dept: Courtroom 1

22  
23  
24 I, Mat Longiaru, declare as follows:

25 1. I am employed by plaintiff Gen-Probe Incorporated as Vice President, Diagnostic  
26 Development. I have been employed by Gen-Probe since February 1991.

27 2. As disclosed in my *Curriculum Vitae* attached hereto as Exhibit 1, I received a B.S.  
28



(Biology) from the City College of New York in 1975, an M.S. (Microbiology) from Long Island University in 1977, and a Ph.D. (Microbiology and Immunology) from Albert Einstein College of Medicine in 1981.

3. As a result of my education and experience, I am familiar with methods of nucleic acid target capture and amplification. I understand methods of non-specific amplification as disclosed in the examples by U.S. Patent No. 5,750,338 ("the '338 patent") and methods of specific amplification, such as Gen-Probe's patented Transcription-Mediated Amplification (TMA) process.

4. I have read the Scientific Background section of the accompanying memorandum in support of summary judgment. The Scientific Background section presents an accurate summary of information about the nucleic acid methods discussed therein.

TMA Uses Sequence-Specific Primers to Achieve Specific Amplification

5. Gen-Probe's HIV-1/HCV Assay ("the Blood Screening Assay") detects small quantities of HIV (human immunodeficiency virus) and HCV (hepatitis C virus) in blood by capturing the viral nucleic acids (i.e., the target nucleic acids) from a sample of blood and amplifying them. The Blood Screening Assay incorporates Gen-Probe's patented TMA technology to *specifically* amplify the captured viral nucleic acids

6. Gen-Probe's Blood Screening Assay achieves specific amplification in part by employing sequence-specific primers, which are designed and made to bind only to specific sequences of interest in the target HIV and HCV nucleic acids. The TMA process will only amplify nucleic acid captured from a sample if the primers find and bind to their respective specific target sequences.

7. One of the two enzymes used in Gen-Probe's Blood Screening Assay is reverse transcriptase ("RT"). Reverse transcriptase is a DNA polymerase that produces a complementary DNA strand copy of a single-stranded RNA or DNA that has a bound primer. In TMA, reverse transcriptase produces complementary DNA from the target nucleic acids (or their complementary strands) only if the sequence-specific primers (described in paragraph 6) first bind to a single strand of RNA or DNA. That is, if the target organism is not present in the sample, the primers



1 will be unable to bind to the captured sequence and the reverse transcriptase will not initiate  
2 synthesis. Therefore the action of the RT enzyme is dependent on the specific primers.

3 TMA Also Uses Specific Promoters and Enzymes to Achieve Specific Amplification

4 8. In addition to its target-specific sequence, some of the TMA primers used in the  
5 Blood Screening Assay contain a "promoter" sequence that allows a specific enzyme, an RNA "T7  
6 polymerase," to *specifically* bind to and produce RNA copies of the target nucleic acids as part of  
7 the TMA amplification process.

8 9. A functional "T7 promoter" is formed in the course of the TMA process if, and only  
9 if, the primer finds and binds to its complementary target sequence in the captured target molecule  
10 so that the target sequence is copied by reverse transcriptase. If the T7 promoter *is* formed as a  
11 result of primer binding to the target sequence, then the T7 RNA polymerase used in Gen-Probe's  
12 Blood Screening Assay will amplify the sequence attached to the T7 promoter sequence. The T7  
13 RNA polymerase does not amplify other sequences present in the sample because they are not  
14 attached to a T7 promoter sequence. Thus, in the Blood Screening Assay, the T7 polymerase  
15 enzyme *specifically* recognizes the T7 promoter sequence, which has been *specifically* attached to  
16 the target sequence by the binding of *specific* primers, and the T7 polymerase *specifically*  
17 amplifies only that sequence.

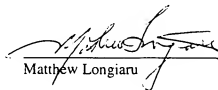
18 10. In the TMA process, one of the primers *specifically* binds to the newly transcribed  
19 RNA and reverse transcriptase makes a new complementary DNA copy of that RNA. The process  
20 repeats in a cyclic fashion, resulting in exponential amplification only of the particular target  
21 sequence of interest as a consequence of the use of sequence-specific primers, specific promoter  
22 sequences, and specific RNA polymerase enzymes. This process safeguards against amplification  
23 of non-target sequences and thus protects against false positive results.

24 11. The TMA method used in the Blood Screening Assay differs substantially from the  
25 non-specific amplification methods disclosed in the '338 patent. All of the methods described in  
26 the examples of the '338 patent *non-specifically* amplify any nucleic acids captured from the  
27 sample, whether those nucleic acids are the intended target or are some other nucleic acid present  
28 in the sample after target capture. Unlike the non-specific amplification methods described in the



1 '338 patent, the TMA process will not amplify non-target sequences that might be retained on the  
2 solid support after the target capture step. The sequence-specific primers, specific promoters, and  
3 specific RNA polymerase enzymes used in TMA are designed to only amplify their intended target  
4 nucleic acids, even if other sequences are present.

5 I hereby declare under penalty of perjury under the laws of the United States of America  
6 that all statements made herein of my own knowledge true and that all statements made on  
7 information and belief are believed to be true. This declaration was executed by me on this 24  
8 day of April, 2001 at San Diego, California.

9  
10   
11 Matthew Longiaru  
12  
13  
14  
15  
16  
17  
18  
19  
20  
21  
22  
23  
24  
25  
26  
27  
28